


# Long Term Care

10 Year Warranty


Gentle Heel Slope and  
Foot Section Channel Cuts  
for Targeted & Effective  
Heel Protection

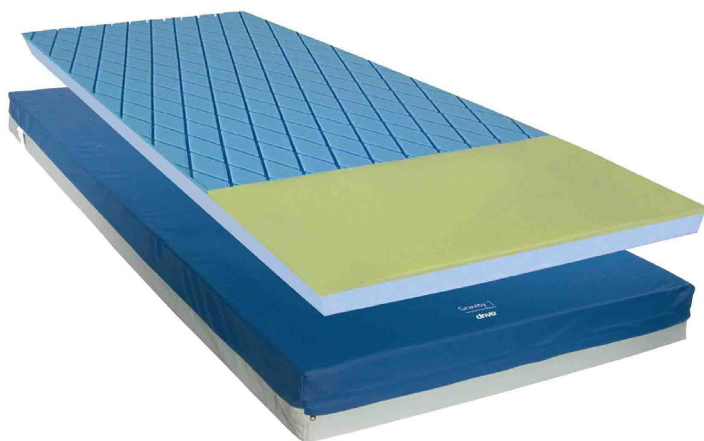

## accessories

15870C	Replacement Cover for 15870	1/ea
15884C	Replacement Cover for 15884	1/ea
15885C	Replacement Cover for 15885	1/ea
15886C	Replacement Cover for 15886	1/ea
15887C	Replacement Cover for 15887	1/ea
15888C	Replacement Cover for 15888	1/ea
158709C	Replacement Cover	

## Heel Slope

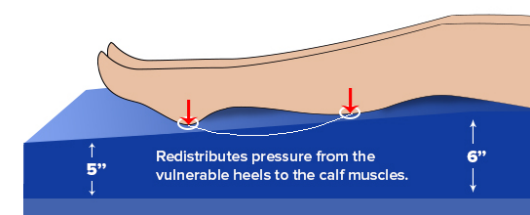

## Gravity 8 Deluxe Long Term Care Pressure Redistribution Mattress

Item #	Raised Side Rails	Product Description	Product Dimensions	Product HCPCS	UOM
15886	No	Multi-Layered/Multi-zoned Foam Mattress	36" (W) X 76" (L) X 6" (H)	E0184	1/ea
15870	No	Multi-Layered/Multi-zoned Foam Mattress	36" (W) X 80" (L) X 6" (H)	E0184	1/ea
15884	No	Multi-Layered/Multi-zoned Foam Mattress	36" (W) X 84" (L) X 6" (H)	E0184	1/ea
15887	Yes	Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	36" (W) X 76" (L) X 6" (H)	E0184	1/ea
15888	Yes	Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	36" (W) X 80" (L) X 6" (H)	E0184	1/ea
15885	Yes	Multi-Layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	36" (W) X 84" (L) X 6" (H)	E0184	1/ea

- Outstanding combination of multi-layered die-cut and memory foam for therapeutic pressure redistribution
- Mattress includes a 6" to 5" heel slope and a "channel cut" center layer of foam under the foot section to target heel ulcer prevention
- Top layer uses unique die-cutting and a true memory foam heel section to optimize zoned pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- Stretch cover is fire retardant, fluid and tear resistant, low shear and easy to clean
- High Density Foam Perimeter provides solid edge support and facilitates safe patient transfers
- Available with or without raised side rails

## specifications

Cover Material	Stretch, Fluid-Resistant Nylon Top   Waterproof Vinyl Bottom
Mattress Layers of Foam	3
Side Rails	Yes
Dedicated Heel Section	Yes
Fire Rating	CFR 16 Part 1633
Weight Capacity	400 lbs.
Warranty	10 Year Limited Foam   1 Year Cover

\*Direct Supply® and Panacea® are registered trademarks of Direct Supply, Inc.